



FaceFX™
ANTI-AGING

Mode d'emploi

Instructions for Use



Treating with Silk'n FaceFX for the First Time

The skin should be, clean, dry and free of any powder, gel or cream.

- Remove Silk'n FaceFX *APPLICATOR* from its cradle in the *BASE UNIT*.
- Press the *POWER ON/OFF SWITCH*. The green *READY/TREATMENT INDICATOR LIGHT (I)* will start blinking, and shortly afterwards will be lit steadily.
- Your FaceFX is ready for treatment.
- Place the *TREATMENT SURFACE* on your skin. The *LED* light will automatically activate.
- Apply light pressure while slowly moving the FaceFX *APPLICATOR* in small circular motions. The *TREATMENT SURFACE* temperature will slowly increase.
- Perform the treatment for 5-7 minutes on each area.
- It is recommended to apply a small amount of a moisturizing cream on the treated skin areas after each session and to massage it gently in circular motions.



Table of Contents

Before You Start	4
What is Silk'n FaceFX?.....	4
Silk'n FaceFX device description.....	4
Package contents.....	5
Safety with Silk'n FaceFX	5
When to avoid using Silk'n FaceFX?.....	5
Reducing the risk of injury.....	6
Possible side effects using FaceFX.....	7
Protect your Silk'n FaceFX.....	8
Skin Treatment the FaceFX Way	9
The benefits of Silk'n FaceFX.....	9
Silk'n FaceFX facial treatment areas.....	9
First Steps with Silk'n FaceFX	9
Silk'n FaceFX setup.....	9
Silk'n FaceFX battery charging.....	9
Plan your Silk'n FaceFX treatment for best results.....	10
Treating with Silk'n FaceFX for the first time.....	11
What to expect when treating with Silk'n FaceFX?.....	12
After treating with Silk'n FaceFX.....	12
Maintenance of Silk'n FaceFX	13
Cleaning Silk'n FaceFX.....	13
Storage.....	13
To protect the environment.....	13
Removing the rechargeable batteries.....	13
Troubleshooting	14
Customer Service	16
Specifications	17
Labeling	18
Français	19

Before You Start

Before using Silk'n FaceFX for the first time, please read this Manual in its entirety. Pay particular attention to sections on device use procedures, device operation, and after-use procedures.

What is Silk'n FaceFX?

Silk'n FaceFX is an FDA cleared hand held light emitting diode (LED) device for facial treatments of periorbital wrinkles, reducing pore size and temporary improvement of the appearance of the skin. The light emitted by Silk'n FaceFX stimulates deep skin tissue and achieves facial rejuvenation.


Typical areas for Silk'n FaceFX treatment are the neck, the cheeks, the forehead, and around the eyes.

Silk'n FaceFX device description

Silk'n FaceFX device consists of a *BASE UNIT*, an *APPLICATOR* and an *ADAPTOR*.

Silk'n FaceFX *APPLICATOR* is a hand held unit that contains *TREATMENT SURFACE* at the *APPLICATOR* tip, *ON/OFF SWITCH* and *INDICATOR PANEL*.

Silk'n FaceFX *BASE UNIT* is used to charge the *APPLICATOR* batteries.

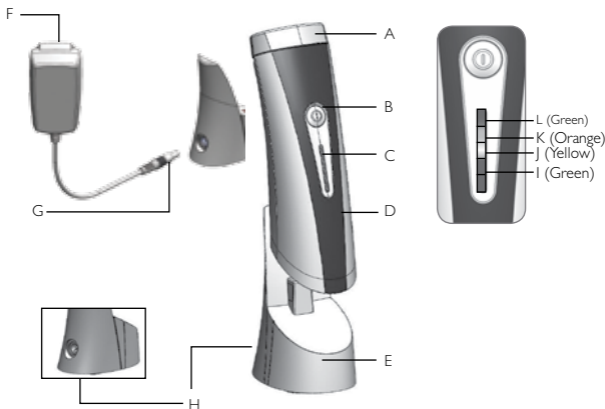

A	Treatment surface	G	Plug
B	On/Off Switch	H	BASE UNIT socket
C	Indicator Panel	I (Green)	Lit-Ready for treatment, Blink-Battery Low
D	APPLICATOR	J (Yellow)	Error (reset the unit and tries again)
E	BASE UNIT	K (Orange)	Battery charging
F	Adaptor	L (Green)	Battery full

Package contents

Upon opening the Silk'n FaceFX package, you will find the following:

- FaceFX *BASE UNIT* and *APPLICATOR*
- Adaptor
- Post treatment cream
- User manual

Safety with Silk'n FaceFX

To protect your skin, a temperature sensor is embedded in the treatment tip of the Silk'n FaceFX. This sensor ensures the skin's surface will not exceed 41°C / 106°F.

When to avoid using Silk'n FaceFX?

Certain conditions may temporarily prevent the use of Silk'n FaceFX. **DO NOT USE Silk'n FaceFX** if any of the following conditions apply to you:

- If you have dark brown or black spots, such as large freckles, birth marks, moles or warts on the area to be treated.
- If you have eczema, psoriasis, lesions, open wounds or active infections, such as cold sores in the area to be treated. Wait for the effected area to heal before using Silk'n FaceFX.
- If you have abnormal skin conditions caused by diabetes or other systemic or metabolic diseases.
- If you are currently or have recently been treated with Alpha-Hydroxi Acids (AHAs), Beta-Hydroxi Acids (BHAs), Retin-A®, topical retinoids or azelaic acid.
- If you have been treated with Accutane® (isotretinoin) within the past 6 months.
- If you have been on a steroid regimen within the past 3 months.
- If you have a history of herpes outbreaks in the area of treatment, unless you have consulted your physician and received preventative before using Silk'n FaceFX

Reducing the risk of injury

- Keep Silk'n FaceFX out of the reach of children.
- Keep Silk'n FaceFX out of the reach of persons with reduced physical, sensory or mental capabilities.



- Keep Silk'n FaceFX away from water!

Do not use or store Silk'n FaceFX in a wet environment (e.g. near a shower or a bathtub)!

Keep the device and its adaptor dry!

- Do not use Silk'n FaceFX if the *ADAPTOR* or the *APPLICATOR* is damaged!
- If the *APPLICATOR* is damaged, make sure it is replaced with an original part from Home Skinovations!
- The *APPLICATOR* contains a dedicated transformer. Never replace the *APPLICATOR* with another plug!

If you are unsure if FaceFX is safe for you to use, please consult with your doctor or dermatologist.

Possible side effects using FaceFX

When used according to the instructions, side effects and complications associated with use of Silk'n FaceFX are rare. However every cosmetic procedure, including those designed for home use, involve some degree of risk. Therefore it is important that you understand and accept the risks and complications that can occur with energy based systems designed for home use.

Silk'n FaceFX treatments are generally very well-tolerated. Most users sense the moderate heat applied by the device during use, with no further discomfort.

• Minor Skin Discomfort

Your skin may become slightly irritated after using Silk'n FaceFX. This is a typical and transient phenomenon. If the irritation does not disappear within three days, please consult a physician or a dermatologist.

• Skin Redness

Your skin may become slightly red right after using FaceFX. This should disappear a few hours after treatment.

Following the treatment it is recommended to apply a small quantity of the Silk'n FaceFX Cream.

Protect your Silk'n FaceFX

- Store the device in a clean area free of dust.
- Do not cover the ventilation slots during use.
- Never shake or drop Silk'n FaceFX or subject it to severe shocks.
- Silk'n FaceFX is sensitive to environmental temperature during use. Do not expose the FaceFX to environmental temperatures lower than 15°C / 59°F or higher than 35°C / 95°F during use.
- Silk'n FaceFX can not tolerate extreme change of temperature. When moved from a cold environment (under 10°C / 50°F) to a warm environment, avoid use for approximately 3 hours.
- Do not expose Silk'n FaceFX to direct sunlight or UV light for long periods of time.
- Silk'n FaceFX is designed for occasional home use (as opposed to exhaustive commercial or professional use). Avoid such exhaustive use.
- Avoid covering Silk'n FaceFX and its adaptor during recharge.

Skin Treatment the FaceFX Way

The benefits of Silk'n FaceFX

The outcome of using Silk'n FaceFX when following instructions correctly will be the reduction of facial fine lines, reducing pore size and improving the appearance of the face. Best results can be expected 3-7 weeks after completing the recommended course of treatment.

Silk'n FaceFX facial treatment areas

Typical areas for Silk'n FaceFX treatment are the cheeks, forehead, neck, and around the eyes.



First Steps with Silk'n FaceFX

Silk'n FaceFX Setup

1. Remove Silk'n FaceFX *BASE UNIT*, *APPLICATOR*, *ADAPTOR* and other components from box.
2. Before using the Silk'n FaceFX for the first time, fully charge the battery. This may take up to 6 hours.

Silk'n FaceFX Battery Charging

- Insert the *APPLICATOR* into its cradle in the FaceFX *BASE UNIT*.

- Plug the *ADAPTOR PLUG* into the Base Unit socket.
- Plug the *ADAPTOR* into an electrical outlet. An orange *BATTERY – CHARGING - INDICATOR LIGHT (K)* will illuminate to indicate that your Silk'n FaceFX is in the process of charging.

Caution!

Never cover Silk'n FaceFX and its adaptor while recharging!

Note!

Silk'n FaceFX can not be used while recharging



- When *BATTERY – CHARGING - INDICATOR LIGHT (K)* turns off and *BATTERY – FULL - INDICATOR LIGHT (L)* illuminates, your Silk'n FaceFX is fully charged and ready to start.

(Recharging will be needed when *BATTERY – CHARGING - INDICATOR LIGHT (K)* illuminates again.)

Plan your Silk'n FaceFX treatment for best results

The efficiency of the results varies from person to person and according to the frequency of use. Some may see results sooner than others.

A single FaceFX session usually takes 15-20 minutes in which you should gently glide the FaceFX applicator slowly over the desired treatment area. During the procedure you should feel gentle heat penetrating deep into the skin.

Initial Silk'n FaceFX treatment course: For the first month (4 weeks), treatment should be performed 3 times a week for 15-20 minutes each time.

Silk'n FaceFX maintenance regime:

Use Silk'n FaceFX once a week for 15-20 minutes.

Treating with Silk'n FaceFX for the First Time

The skin should be, clean, dry and free of any powder, gel or cream.

- Remove Silk'n FaceFX *APPLICATOR* from its cradle in the *BASE UNIT*.
- Press the *POWER ON/OFF SWITCH*. The green *READY/TREATMENT INDICATOR LIGHT (I)* will start blinking, and shortly afterwards will be lit steadily.

Your FaceFX is ready for treatment.

- Place the *TREATMENT SURFACE* on your skin. The *LED* light will automatically activate.
- Apply light pressure while slowly moving the FaceFX *APPLICATOR* in small circular motions. The *TREATMENT SURFACE* temperature will slowly increase.
- Perform the treatment for 5-7 minutes on each area.
- It is recommended to apply a small amount of a moisturizing cream on the treated skin areas after each session and to massage it gently in circular motions.



What to expect when treating with Silk'n FaceFX?

During a Silk'n FaceFX session it is normal to experience and feel:

- **A Sensation of Warmth** – During the Silk'n FaceFX session you may feel slight warmth. This is normal and should disappear up to an hour after the session.
- **Some Mild Red Color** – After each Silk'n FaceFX session you may notice a slight reddening of the skin (mild erythema) over the treated area. This is normal and should fade within an hour or two.

Best results can be expected 3-7 weeks after completing the recommended course of treatment.

After treating with Silk'n FaceFX

- When FaceFX session has been completed turn FaceFX off by pressing the *POWER ON/OFF SWITCH*.
- After each session it is recommended that you clean your FaceFX device, especially the *APPLICATOR TIP* (See Section: "Cleaning FaceFX").
- Place the *APPLICATOR* in its *BASE UNIT* cradle to recharge the batteries (See Section: "Silk'n FaceFX Battery Charging").

Cleaning the Silk'n FaceFX

After each session, it is recommended to clean your Silk'n FaceFX device:

- Turn Silk'n FaceFX off by pressing the *POWER ON/OFF SWITCH*.
- Unplug Silk'n FaceFX from the electrical outlet
- Allow cooling down of its *TREATMENT SURFACE*.
- Use a dry, clean cloth and a specially formulated cleaner for electronic equipment to gently wipe the device and especially its *TREATMENT SURFACE*.



- Never immerse Silk'n FaceFX or any of its parts in water!
- Never clean the FaceFX or any of its parts under the tap or in the dishwasher.

Storage

- Store the clean device and all its parts in a dry place, at a temperature of between -10°C to 60°C / 14°F to 140°F

To protect the environment

At the end of its life, do not discard Silk'n FaceFX with the normal household waste. Dispose of it at an official recycling collection point.



- Remove batteries and dispose at an official collection point for batteries. Please call your local recycling depot for locations and proper recycling instructions.

Removing the rechargeable batteries

Remove the rechargeable batteries only if they are completely empty and cannot be recharged.



Undo the screw at the bottom panel of the *APPLICATOR*.



Lift and remove the top panel of the *APPLICATOR*.



Push the three snaps and pull out the electronics card.



Cut the tow wires that connect the rechargeable batteries pack to the *APPLICATOR*.



Cut the tow wires that connect the electronics card to the lower panel of the *APPLICATOR*.

Troubleshooting

This chapter summarizes the most common problems you could encounter. If you are unable to solve the problem with the information below, contact your local Silk'n Customer Service Center for further support.

Problem	Possible Cause	Solution
Silk'n FaceFX does not work	The rechargeable batteries are empty	Charge the <i>APPLICATOR</i>
	The <i>APPLICATOR</i> is broken	Contact your local Silk'n customer Service center
Silk'n FaceFX suddenly switches off	The rechargeable batteries are empty	Charge the <i>APPLICATOR</i>
Device temperature rises during recharge	This is normal	No action is required

BATTERY – CHARGING - INDICATOR LIGHT (K Orange) does not turn on when connecting the PLUG into the BASE UNIT SOCKET	The ADAPTOR is not plugged in the electrical outlet	Insert the <i>PLUG</i> into the <i>BASE UNIT SOCKET</i> Connect another appliance to the wall socket to check if the electrical outlet is live If the <i>APPLICATOR</i> still does not charge, contact your local Silk'n Customer Service Center
	The ADAPTOR is not inserted properly into the electrical outlet and/or the <i>PLUG</i> is not inserted properly into the <i>BASE UNIT SOCKET</i>	Make sure the <i>APPLICATOR</i> is inserted properly into the electrical outlet and/or <i>PLUG</i> is inserted properly into the <i>BASE UNIT SOCKET</i>
	The <i>APPLICATOR</i> is broken	Contact your local Silk'n Customer Service Center
APPLICATOR temperature rises during use	This is normal	No action is required
READY/TREATMENT MODE-INDICATOR (I Green) does not lit when pressing the ON/OFF SWITCH	Batteries are empty	Charge the applicator
	Batteries are full (I Green) indicator lit when applicator is placed in base unit charger	Press the <i>ON/OFF SWITCH</i> twice to switch the applicator off and on again. If <i>READY/TREATMENT MODE-INDICATOR (I Green)</i> does not lit again, contact your local Silk'n Customer Service Center
The Yellow SYSTEM ERROR - INDICATOR LIGHT (J Yellow) is lit.	Malfunction of the <i>APPLICATOR</i>	Press the <i>ON/OFF SWITCH</i> twice to switch the <i>Applicator</i> off and on again. If <i>SYSTEM ERROR - INDICATOR LIGHT (J Yellow)</i> illuminates again, contact your local Silk'n customer Service center
The Yellow SYSTEM ERROR - INDICATOR LIGHT (J Yellow) is blinking	The rechargeable batteries are empty	Charge the <i>ADAPTOR</i>
The Orange BATTERY – CHARGING and Green BATTERY – FULL - INDICATOR LIGHTS are illuminated alternately	Malfunction of the <i>ADAPTOR</i>	Disconnect the <i>ADAPTOR</i> from the electrical outlet, press the <i>ON/OFF SWITCH</i> twice to switch the <i>APPLICATOR</i> on and off again. Connect the <i>Adaptor</i> to the wall socket, If the problem persists, contact your local Silk'n Customer Service Center
	Malfunction of the <i>APPLICATOR</i>	Press the <i>ON/OFF SWITCH</i> twice to switch the <i>Applicator</i> off and on again. If the problem persists, contact your local Silk'n Customer Service Center

Customer Service

For USA and Canada customers:

For more information about Silk'n™ products please visit www.silkn.com.

If your device is broken, damaged, in need of repair, or for any other user assistance, please contact Silk'n™ Customer Service: **1-877-DO-SILKN / 1-877-367-4556, contact@silkn.com**

International Customers:

For more information about Silk'n products or to find a list of International Silk'n websites please visit www.silkn.eu or www.silkn.com.

If your device is broken, damaged, in need of repair, or for any other user assistance, please contact your local Silk'n™ Customer Service or write to: info@silkn.com.

Manufacturer:


Home Skinovations Ltd.

Apolo Building, Shaar Yokneam POB 533


Yokneam 20692





ISRAEL

Specifications

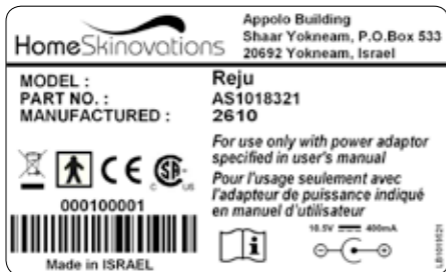Rated voltage (Adaptor)	100V - 240V
	Model: UE08WCP-105040SPA
Rated frequency (Adaptor)	50Hz – 60Hz
Rated current (Adaptor)	400 mA
Applied part	 Type BF
Protecting rating	IP20 (EN60529)
Storage condition	Temperature: -10°C to 60°C / 14°F to 140°F Relative humidity: 5% to 95%
Operating condition	Temperature: 15°C to 35°C / 59°F to 95°F Relative humidity: 10% to 80% Atmospheric pressure: 700 to 1060 hPa
Light source	Light Emitting Diodes Wavelengths peak emission wavelength: 645nm Intensity power: 250mW ±35%
Lithium ion battery	2 x 3.7V, 750±50mAh
Max. charging time	6 hours


Labeling



Degree of protection against electric shock type BF applied part	
WEEE - Waste Electrical and Electronic Equipment	
CSA Mark	
CE Mark	

Degree of protection against ingress water: ordinary

This device is not suitable for use in the presence of flammable anesthetic mixture with air or with Oxygen or Nitrous Oxide.

The symbol  on the label affixed to this device means “Attention, consult accompanying documents”.



Apolo Building, Shaar Yokneam POB 533
Yokneam 20692, ISRAEL

tel. +972-4-9097470, fax +972-4-9097471
info@silkn.com, contact@silkn.com

www.silkn.com



PB101950A